

DA0970

Specimen Preservation Reagent



Intended use

This kit is used for the preservation of tissue and cytopathological analysis specimen.

Features



Virus inactivation



Prevent RNA degradation




Transportion and storage at RT

Principle

Guanidinium is the main component of the Specimen Preservation Reagent. It can lyse and denature virus and microorganisms in the cells, meanwhile it inhibits/denature the RNase, thereby preserving the RNA in the specimen.

Procedure

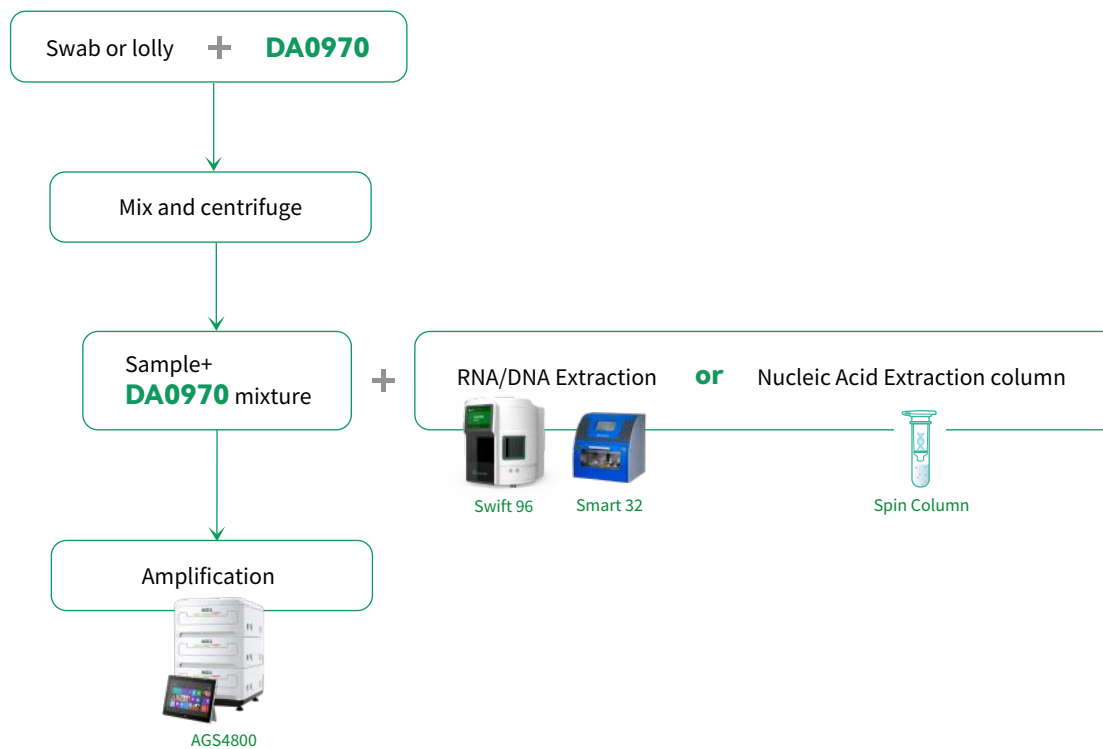

Order Information

Description	Cat. No	Specification(Tests/Kit)
Specimen Preservation Reagent	DA0970	1.5mL/tube, 100tubes/kit;
		3mL/tube, 100tubes/kit;
		5mL/tube, 100tubes/kit;
		10mL/tube, 50tubes/kit;
		500mL/bottle, 1 bottle/kit

Daan Gene Co., Ltd.

CHINA (Headquarters)

Add: No.19, Xiangshan Road, Science Park, High & New Technology Development District, Guangzhou, Guangdong, P.R. China

Tel.: +86-20-32068126

Fax: +86-20-32068352

Email: marketing@daangene.com

CANADA office

Add.: #413-4538 Kingsway, Burnaby, BC V5H 4T9, Canada

Tel.: +1604 451 7588

Fax: +604451.7587

Email: gbai@daangene.net

Website: <https://en.daangene.com/>

